



Shanghai Tianfan Pharmaceutical Machinery Factory Events

- *1978—Shanghai Tianma Pharmaceutical Components Factory (collectively-owned enterprise)
- *1980—Shanghai Tianma Pharmaceutical Machinery Factory (collectively-owned enterprise);
- *1999--Restructuring of enterprise;
- *2000—The senior engineers established Shanghai Tianfan Parmacuetical Machinery;
- *2005—ISO-9001;
- *2007—EU CE;
- *2015—Moved to new factory, building area :15300sqm;
- *2016—Mold processing center established.



WELCOME TO US




michael@stctp.com whatsapp+8618017027102


Factory brief video



CNC center, must for high quality & accuracy

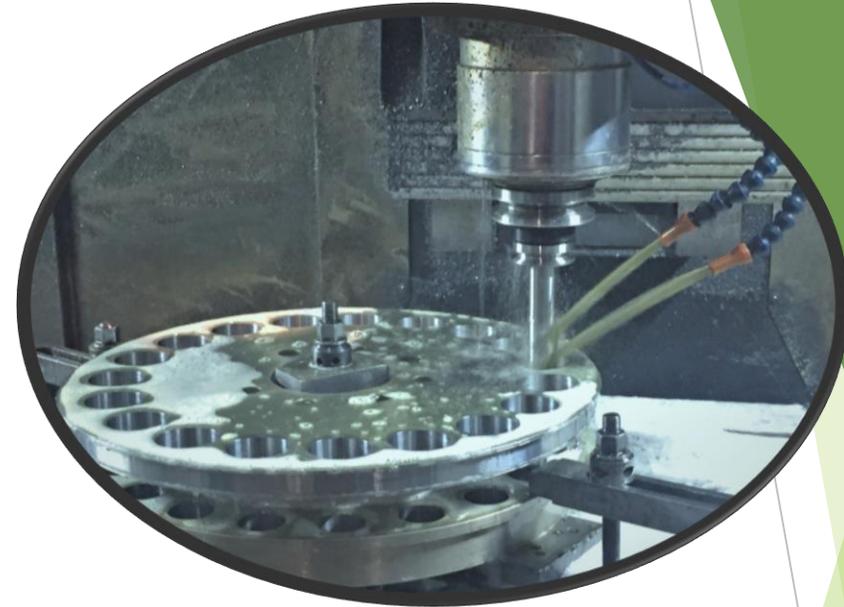

Molds workshop providing customizing service

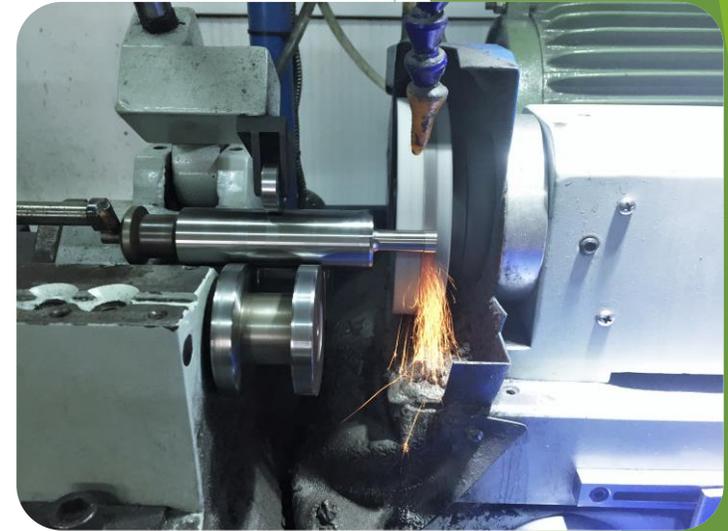

Finished machines waiting for QC test



Assembling workshop 1



Assembling workshop 2



Packing up area



ZPT420 Series upgraded model





ZPT226 Series upgraded model





GZPK370 Series high-speed tablet press





GZPK370i Series Semi-automatic tablet press





GZPK370i Series automatic tablet press





What is ZPT model?

ZPT models are based on normal ZP models but having some functions of automatic rotary tablet press!

Main features of ZPT are as below:

ZPT model has dust-prevent ring, oil catching rings

1. Dust-prevent ring:

It can prevent powder going inside the upper mold holder and oil flowing down the punch rods, reduce noise and abrasion.


2. Oil catching ring:

It can catch the lubrication oil if any leaking happens.

3. Side fixed bolt:

Part of the bolt is in the side of the middle die, which can prevent the bolt coming off effectively.

ZPT model

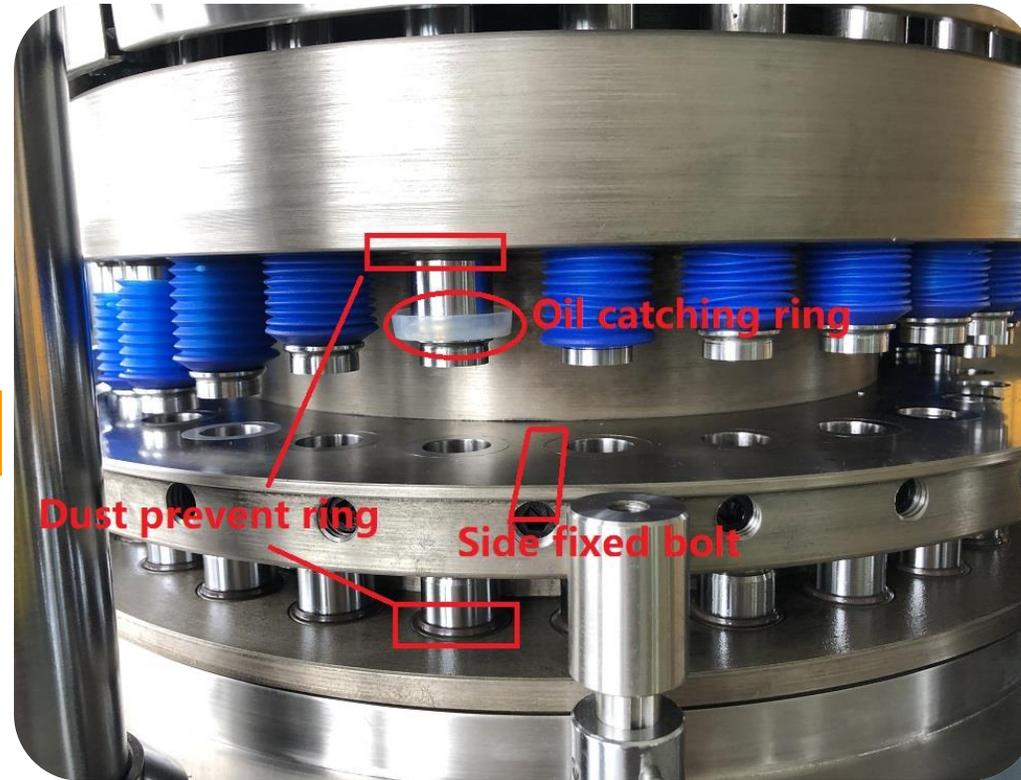

ZP model



ZPT lubrication pump and used oil collection system

auto-lubrication pump

used oil collection system



部分厂家虽安装自动油
泵，但并未增设废油回收
系统，导致废油长期滞留
转盘内部，吸附粉尘后严
重影响设备寿命!!!

Used lubrication oil container
废油回收盒



ZPT used oil guiding slots and holes


ZPT turret

Guiding slots:

The guiding slots and holes can help the used oil flow down into the waste oil container. The dust will pile up inside the turret if without these slots, which surely will cause trouble to the turret in the long run.



ZP turret

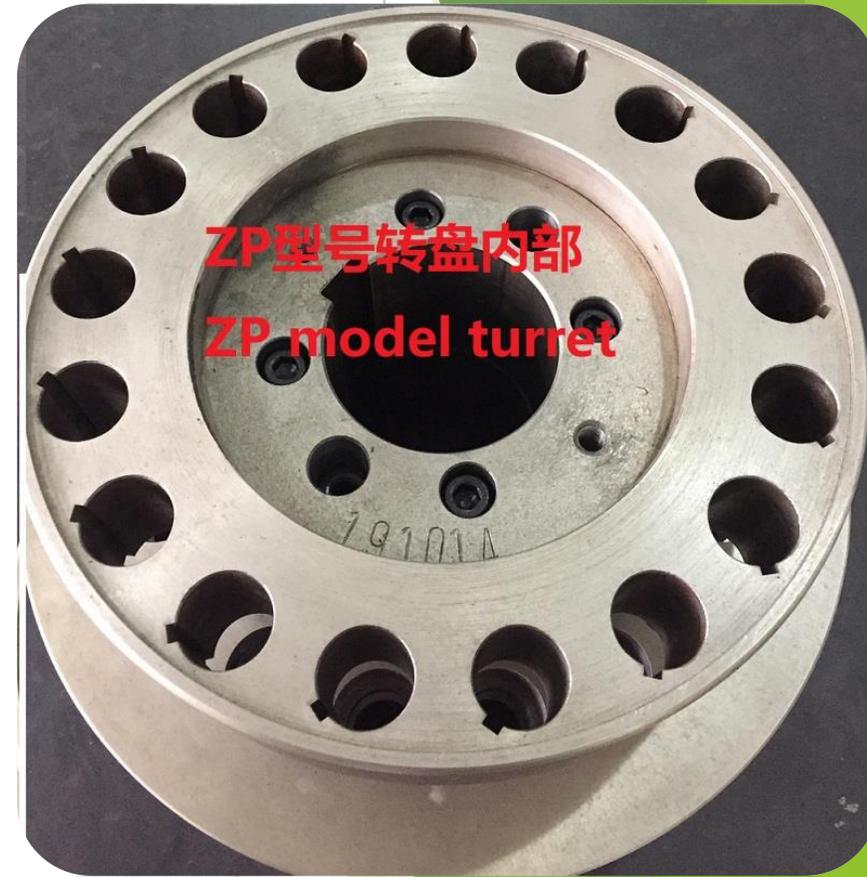


ZPT型号自动润滑系统***油刷

Lubrication brush for upper punches



Lubrication brush for down punches

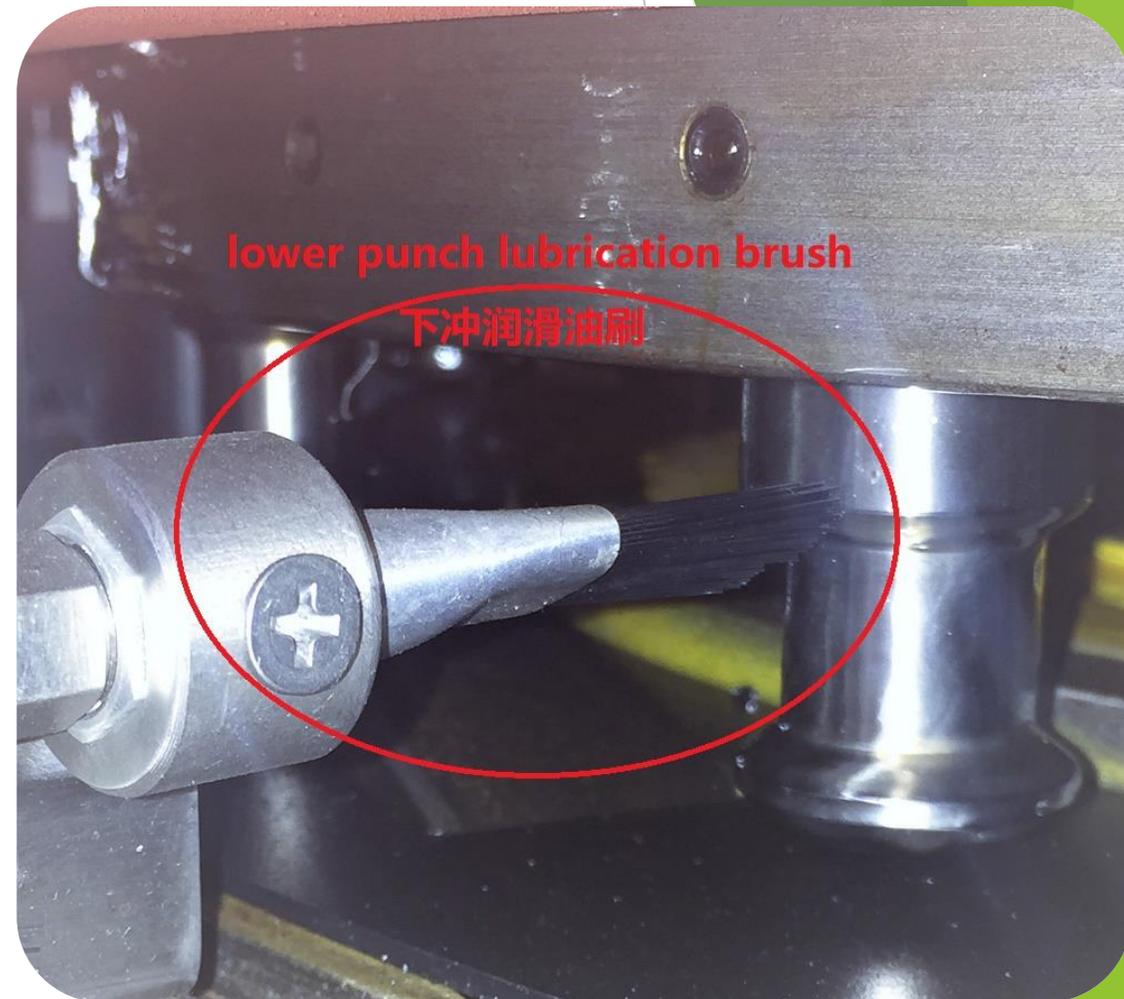

ZPT molds longer than ZP punches



ZPT

ZP

IPT



Thank you for your time!



®

天凡药机
泰博雷特



Sweep and add me WeChat

Michael Sales Manager

M: +86 18017027102

Shanghai Tianfan Pharmaceutical Machinery Factory
Shanghai Tablet Press Mechanical Components Co.,Ltd.

Add: No. 159 Zhu Jia Bang Rd, Xiaokunshan
Town, Songjiang District, Shanghai

Tel: 86-21-57668183

Fax: 86-21-57668600

Q Q: 2850665073 PC: 201614

E-mail: michael@stctp.com

Service hotline: 400880-5118

Specialized in Powder forming Solutions

www.tablet-press.net

michael@stctp.com
[whatsapp+8618017027102](https://www.whatsapp.com/channel/00299a618017027102)