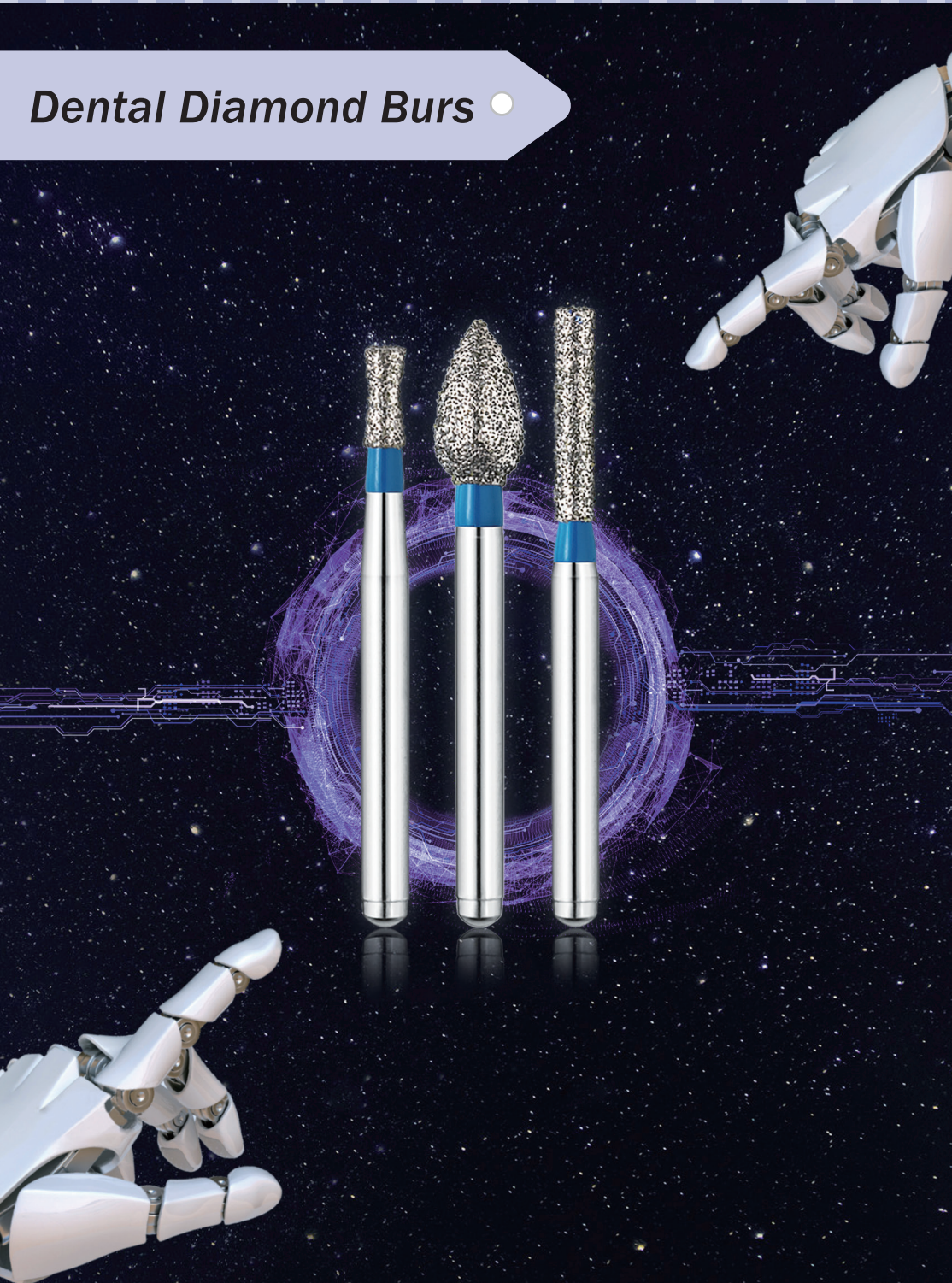
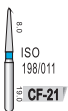
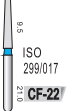
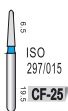
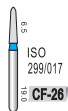
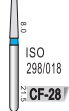
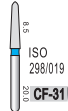
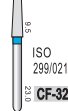



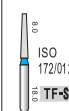
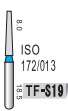
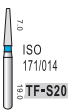


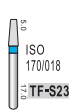
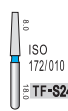

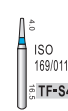
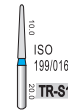

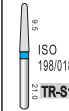
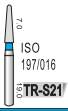
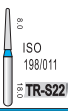
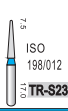
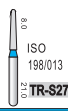
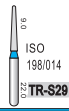

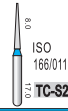

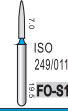
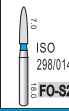
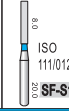
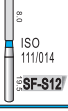
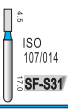
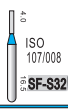

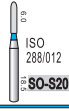

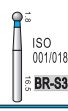

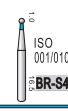














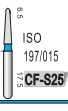





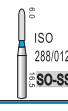

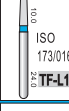
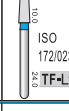
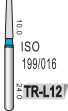
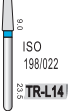
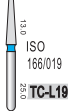

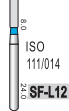
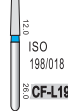
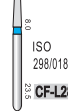

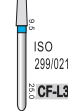
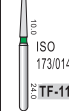

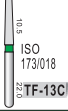
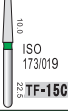
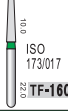
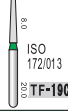

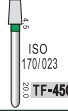
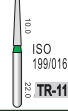


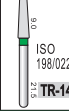
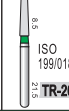
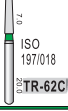
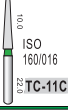
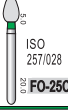
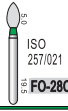




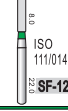
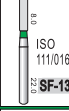
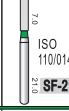
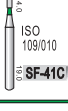
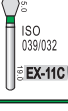







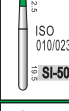
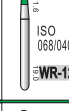



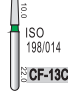







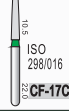


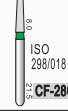
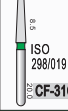
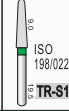
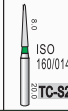
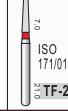
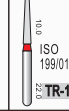
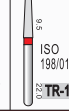
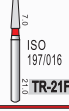
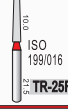
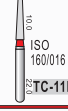
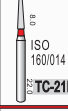

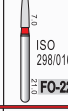
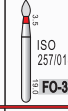
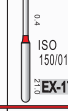







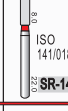
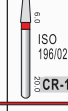

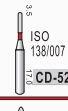

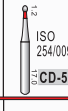
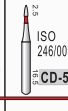
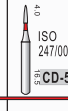
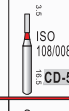

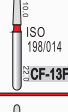

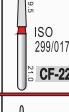
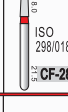
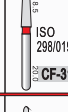
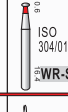
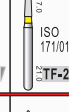




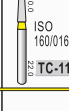
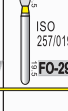
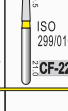
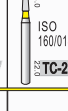

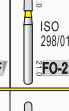










Dental Diamond Burs

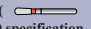


ISO 173/014 TF-11	ISO 173/016 TF-12	ISO 173/018 TF-13	ISO 172/023 TF-14	ISO 173/019 TF-15	ISO 173/017 TF-16	ISO 173/018 TF-17	ISO 172/012 TF-18	ISO 172/013 TF-19	ISO 171/014 TF-20	ISO 171/016 TF-21
ISO 170/021 TF-22	ISO 170/018 TF-23	ISO 198/012 TF-24	ISO 170/016 TF-31	ISO 169/011 TF-41	ISO 170/012 TF-42	ISO 170/014 TF-43	ISO 199/016 TF-11	ISO 199/016 TF-12	ISO 198/018 TR-13	ISO 198/022 TR-14
ISO 200/022 TR-15	ISO 198/018 TR-16	ISO 200/024 TR-19	ISO 197/014 TR-20	ISO 197/016 TR-21	ISO 198/011 TR-22	ISO 198/012 TR-23	ISO 197/018 TR-24	ISO 198/016 TR-25	ISO 199/018 TR-26	ISO 198/013 TR-27
ISO 198/014 TR-29	ISO 299/016 MS-13	ISO 299/019 MS-14	ISO 166/011 TC-10	ISO 166/016 TC-11	ISO 166/011 TC-20	ISO 160/014 TC-21	ISO 160/012 TC-26	ISO 299/012 FO-11	ISO 250/014 FO-12	ISO 250/017 FO-13
ISO 251/023 FO-14	ISO 249/011 FO-15	ISO 249/013 FO-16	ISO 298/014 FO-21	ISO 298/016 FO-22	ISO 249/016 FO-23	ISO 257/028 FO-25	ISO 257/032 FO-27	ISO 257/019 FO-29	ISO 257/018 FO-30	ISO 257/018 FO-32
ISO 540/009 SF-38	ISO 111/012 SF-11	ISO 111/014 SF-12	ISO 111/016 SF-13	ISO 110/014 SF-21	ISO 111/012 SF-22	ISO 109/013 SF-31	ISO 107/008 SF-32	ISO 109/010 SF-41	ISO 131/014 SO-15	ISO 288/012 SO-20
ISO 289/014 SO-21	ISO 130/012 SO-24	ISO 141/010 SR-10	ISO 141/012 SR-11	ISO 141/014 SR-12	ISO 141/016 SR-13	ISO 140/012 SR-18	ISO 140/014 SR-19	ISO 553/019 RS-11	ISO 545/018 RS-21	ISO 544/018 RS-31
ISO 039/032 EX-11	ISO 039/035 EX-12	ISO 237/018 EX-20	ISO 237/012 EX-21	ISO 002/018 EX-23	ISO 220/017 EX-24	ISO 237/032 EX-26	ISO 655/016 EX-27	ISO 655/018 EX-28	ISO 655/021 EX-29	ISO 234/014 EX-31
ISO 237/010 EX-41	ISO 001/023 BR-30	ISO 001/018 BR-31	ISO 001/016 BR-40	ISO 001/014 BR-41	ISO 001/010 BR-45	ISO 001/012 BR-46	ISO 001/008 BR-49	ISO 002/018 BC-31	ISO 002/014 BC-42	ISO 002/012 BC-43
ISO 002/010 BC-44	ISO 002/012 BC-45	ISO 019/010 DI-40	ISO 019/010 DI-41	ISO 019/014 DI-42	ISO 019/010 SI-45	ISO 010/012 SI-46	ISO 010/014 SI-47	ISO 010/016 SI-48	ISO 010/023 SI-50	ISO 068/042 WR-11
ISO 068/040 WR-12	ISO 068/042 WR-13	ISO 068/042 WR-18	ISO 903/043 WR-19	ISO 903/010 CF-11	ISO 198/013 CF-12	ISO 198/014 CF-13	ISO 198/015 CF-14	ISO 198/017 CF-15	ISO 298/016 CF-17	ISO 198/019 CF-18

Color Coding

Color code indication is provided at the end of working part (), according to the size of diamond grit conforming to the I.S.O. specification.

Classification

UADENTAL codina system reflects shapes of burs. The contour and sizes could be imagined by UADENTAL nomenclature since I.S.O. numbers are also indicated.

- C/Coarse
- F/Fine
- EF/Extra Fine
- S/Short Shank
- SS/Super Short Shank
- L/Extra length

Type

- TF/Taper Flat End
- TR/Taper Round End
- TC/Taper Conical End
- FO/Flame Ogival End
- SF/Straight Flat End
- SO/Straight Ogival End
- SR/Straight Round End
- RS/Rounded Shoulder
- EX/Special(Extra)Shape
- BR/Ball Round Type
- BC/Ball Collar Type
- WR/Wheel Round Edge
- SI/Single Inverted Cone
- CD/Children's Dentistry Dia-Burs
- DI/Double Inverted Cone
- CR/CR Inlay Preparation Dia-Burs
- MS/Marginal Shaper
- CF/Chamfur

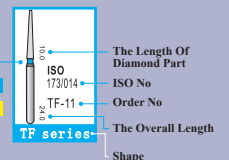
High Quality Control

100 percent inspection with laser measurement machine, cutting test and strength test for each lot are performed. All products comply with our own standards which are far more strict than I.S.O. standards.

BEFORE USE UADENTAL DIA-BURS, PLEASE SEE THE INSTRUCTION AS INDICATED BELOW.

- For dental use only.
- Non-sterile.
- Insert the shank into handpiece according to the instruction of the handpiece manufacturer. Insertion should be the full length.
- Before use, run the instruments to make sure there is no descentering. If the instruments become bent, be sure to remove them from any further use.
- Please always use spray. A plentiful water supply is necessary to maintain long-lived performance and also to prevent loading and piliptis.
- Handle with care not to put excessive load. Pressure should be 20 to 50p for all instruments with diameter up to 0.23.
- Use below max. allowable r.p.m. of each instrument, see chart above, to avoid a possibility of being injured due to breakage.
- Avoid levering, titing, twisting and jamming motion. Do not work with cutting.
- Wear a safety glasses and/or facemask to protect from being injured.
- After using, follow appropriate treatment as the medical waste.
- Practicable for using auto-clave, chemical vaporized-clave, dry heat and E.O.G. sterilization.

How To See This Chart



Color Coding

- Coarse /125-150 μ
- Standard /106-125 μ
- Fine / 5.3 - 6.3 μ
- ExtraFine /20-30 μ